

The FAIR Premium
Implant-System

ICX

ICX-ACTIVE

...the ICX-Implant for immediate restoration



FOR IMMEDIATE RESTORATION:

- ✓ bone-condensing
- ✓ Possibility of active implant orientation
- ✓ optional resection



ICX-ACTIVE

...the ICX-Implant for immediate restoration

A concept for special needs of dentists and clinicians

After losing teeth, patients want immediate and rapidly fixed teeth at affordable costs. Conventional treatments that often involve healing phases that last months on end are often discarded. An immediate restoration with immediate implantation is what is needed.

The traditional two-stage implant protocol with delayed restoration has a long treatment duration, in which patients had to wear a removable appliance and experienced significant discomfort during the recovery period.

Nowadays, immediate implant placement into an extraction site followed by immediate restoration of a dental implant can shorten the dental rehabilitation time and preserve the patients' esthetic appearance at all stages of treatment. medentis now expands the successful ICX implant concept with a new next-generation implant for immediate restoration.

ICX-ACTIVE offers dentists a high degree of flexibility regarding insertion and prosthetic options.


The ICX-ACTIVE implant was designed in order to meet increasing demands of dentists.

- High primary stability even in inconvenient bone conditions
- Bone-condensing characteristics
- active implant orientation for an optimal final implant insertion
- integrated platform switching
- prosthetic connection with double features

coronal ICX-design

variable double thread

ICX-Implant body design

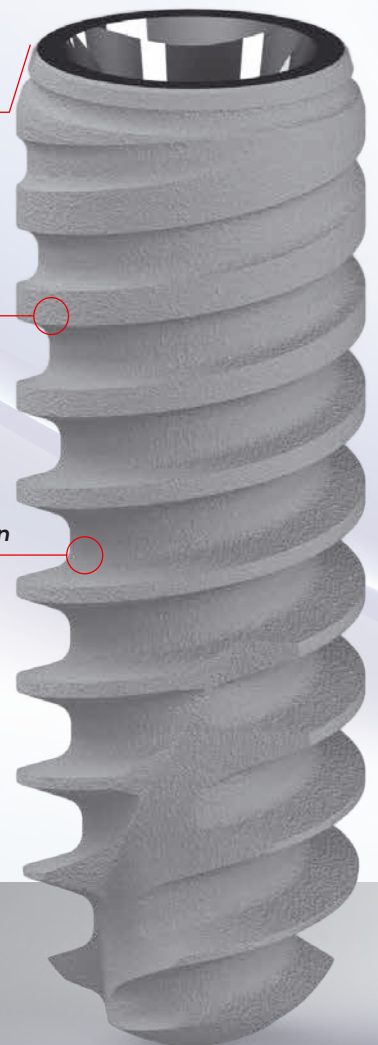


**The FAIR Premium
Implant-System** **ICX**

Service-Tel.: +49 (0)2641 9110-0 · www.medentis.de

An innovative design, which enables immediate restoration with immediate implantation adapted to patients needs

The ICX-Active implant has a tapered implant design with sharp, widely spaced, gradually expanding, double-helical threads that have a steeper thread pitch than most implants.

The implant design gives a high level of primary stability and makes it an excellent choice for immediate extraction, immediate placement, and immediate load cases,

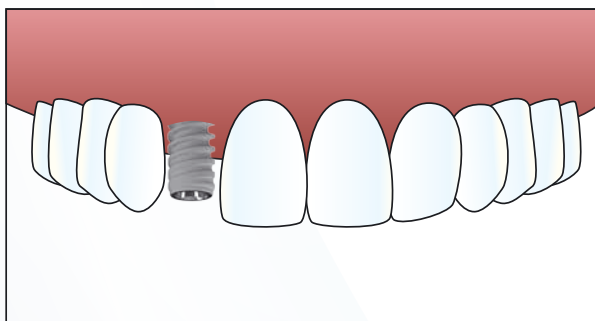


as well as in cases with poor or questionable quality of bone. In addition, these implant design features allow for the redirection of the implant during placement while maintaining high levels of stability.

The platform allows for either "platform switching" or abutment fabrication directly on the platform. These design features will convince users.

The ICX-Active design features are developed to minimize crestal bone loss and maximize the maintenance of gingival architecture. This implant type offers clinicians an excellent solution for restorations in the esthetic zone.

The advanced body design allows clinicians to adjust the implant position for optimal restorative orientation. It is indicated for all positions, but is especially effective in regions of diminished bone quality or quantity, due to its bone-condensing capability,

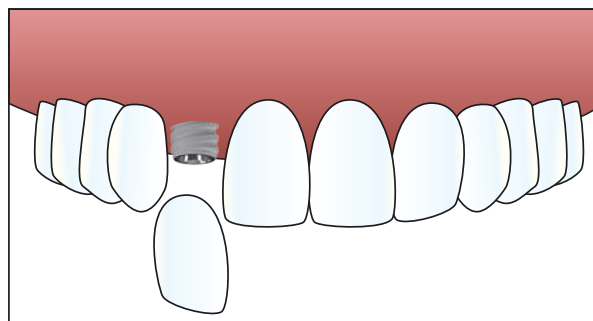

capacity for shorter drilling protocols, and ability to change direction or insertion. It is a viable solution for difficult treatment situations, often providing quicker treatments, shorter healing times, and less invasive procedures. ICX-Active implants are compatible with a wide range of ICX prosthetics and abutments (www.medentis.de).

CLINICAL BENEFITS:

- ✓ rapid insertion ✓ osteotom-effect ✓ increased insertion torque
- ✓ High primary stability even in inconvenient bone conditions
- ✓ Bone-condensing features ✓ integrated platform switching
- ✓ active implant orientation for an optimal final implant insertion
- ✓ prosthetic connection with double features



ICX-ACTIVE

...the ICX-Implant for immediate restoration

X	Article-No.		Number	Price* (Euro)
		ICX-ACTIVE-Implants		
	ACT-001-425100	ICX-ACTIVE Ø 4.2 L 10mm (only available in 6er Packages/Piece € 59,-)		354,-
	ACT-001-425125	ICX-ACTIVE Ø 4.2 L 12.5mm (only available in 6er Packages/Piece € 59,-)		354,-
	ACT-001-425150	ICX-ACTIVE Ø 4.2 L 15mm (only available in 6er Packages/Piece € 59,-)		354,-

Name

Customer-Nr.

Date, Sign

*All Prices in Euro including VAT until Recall (Date 01.12.16)

Praxis Stamp